



ALMA MATER STUDIORUM  
UNIVERSITÀ DI BOLOGNA

ARCHIVIO ISTITUZIONALE  
DELLA RICERCA

Alma Mater Studiorum Università di Bologna  
Archivio istituzionale della ricerca

Design and in-house validation of a portable system for the determination of free acidity in virgin olive oil

This is the final peer-reviewed author's accepted manuscript (postprint) of the following publication:

*Published Version:*

Design and in-house validation of a portable system for the determination of free acidity in virgin olive oil / Grossi M.; Palagano R.; Bendini A.; Ricco B.; Servili M.; Garcia-Gonzalez D.L.; Gallina Toschi T.. - In: FOOD CONTROL. - ISSN 0956-7135. - ELETTRONICO. - 104:(2019), pp. 208-216. [10.1016/j.foodcont.2019.04.019]

*Availability:*

This version is available at: <https://hdl.handle.net/11585/700389> since: 2019-09-25

*Published:*

DOI: <http://doi.org/10.1016/j.foodcont.2019.04.019>

*Terms of use:*

Some rights reserved. The terms and conditions for the reuse of this version of the manuscript are specified in the publishing policy. For all terms of use and more information see the publisher's website.

This item was downloaded from IRIS Università di Bologna (<https://cris.unibo.it/>).  
When citing, please refer to the published version.

(Article begins on next page)

This is the final peer-reviewed accepted manuscript of:

Design and in-house validation of a portable system for the determination of free acidity in virgin olive oil. Grossi M.; Palagano R.; Bendini A.; Ricco B.; Servili M.; Garcia-Gonzalez D. L.; Gallina Toschi T.

FOOD CONTROL, Volume: 104, pages: 208-216

The final published version is available online at:

<https://doi.org/10.1016/j.foodcont.2019.04.019>

© 2019. This manuscript version is made available under the Creative Commons Attribution-NonCommercial-NoDerivs (CC BY-NC-ND) 4.0 International License (<http://creativecommons.org/licenses/by-nc-nd/4.0/>)

1                   **Design and in-house validation of a portable system for the determination**  
2                                   **of free acidity in virgin olive oil**

3  
4                   Marco Grossi<sup>a</sup>, Rosa Palagano<sup>b,\*</sup>, Alessandra Bendini<sup>b</sup>, Bruno Riccò<sup>a</sup>,  
5                   Maurizio Servili<sup>c</sup>, Diego Luis García-González<sup>d</sup>, Tullia Gallina Toschi<sup>b</sup>

6  
7                   <sup>a</sup> Department of Electrical Energy and Information Engineering “Guglielmo Marconi” (DEI),  
8                   University of Bologna, Bologna, 40136, Italy

9                   <sup>b</sup> Department of Agricultural and Food Sciences (DISTAL), University of Bologna, Cesena, 47521,  
10                   Italy

11                   <sup>c</sup> Department of Agricultural, Food and Environmental Sciences, University of Perugia, Perugia,  
12                   06121, Italy

13                   <sup>d</sup> Instituto de la Grasa (CSIC), Sevilla, 41013, Spain

14  
15                   \* Corresponding author: [rosa.palagano@unibo.it](mailto:rosa.palagano@unibo.it) (Rosa Palagano), Department of Agricultural and  
16                   Food Sciences (DISTAL) - University of Bologna, Piazza Goidanich 60, Cesena (FC), 47521, Italy,  
17                   Tel +39 0547 338121.

18 **ABSTRACT**

19 Nutritional and healthy values are well known properties of virgin olive oil (VOO). The product  
20 quality, in terms of belonging to a specific quality grade (extra virgin, virgin, lampante), is defined  
21 by a set of chemical-physical and sensory measurements. According to the official regulation of the  
22 European Union (EU Reg. 1348/2013) the free acidity is the first parameter that has to be  
23 determined by analysts; it gives information about the quality of the olives used to produce the  
24 VOO as well as the hydrolytic state of VOO just produced and stored. The official procedure is  
25 based on an acid-base titration that needs to be carried out in a chemical laboratory.

26 In this paper a portable battery-operated electronic system to measure olive oil free acidity is  
27 presented: the system can be used for quick “in situ” tests in a production environment (olive oil  
28 mills or packaging centers) by people without particular training. The working principle of the  
29 system is based on the creation of an emulsion between oil and a hydroalcoholic solution: the free  
30 acidity is estimated on the value of the emulsion electrical conductance.

31 **The proposed** system has been calibrated and in-house validated showing good results in terms of  
32 limit of detection and quantification, precision and accuracy. **Moreover, a good correlation ( $R^2_{adj} =$**   
33 **0.97) with free acidity data obtained applying the official method on 30 olive oil samples belonging**  
34 **to different commercial categories (extra virgin, virgin and lampante olive oil) has been evidenced.**

35  
36 **Keywords:** portable system; free acidity; virgin olive oil; impedance spectroscopy; in-house  
37 validation; electrical conductance.

38  
39 **Abbreviations**

40 **EIS:** electrical impedance spectroscopy, **EVOO:** extra virgin olive oil, **FA:** free acidity, **LOD:**  
41 limit of detection, **LOO:** lampante olive oil, **LOQ:** limit of quantification, **RSD:** relative standard  
42 deviation, **TAGs:** triglycerides, **VOO:** virgin olive oil.

## 43 1. Introduction

44 Virgin olive oil (VOO) is obtained from olives (the fruits of *Olea europaea L.*) applying only a  
45 mechanical-physical extraction process and represents a product highly appreciated for its  
46 beneficial effects on human health, mainly due to a high content of oleic acid and minor  
47 components such as phytosterols, carotenoids, tocopherols and hydrophilic phenols (Bendini et al.,  
48 2007). The European Commission Implementing Regulation 1348/2013 (EU Reg. No 1348/2013)  
49 defines a decision tree for verifying whether a VOO is consistent with the category declared, and  
50 the quality criteria that have to be checked by analysts are: free acidity (FA), peroxide value,  
51 specific extinctions in UV, sensory characteristics and ethyl esters of fatty acids. The first quality  
52 parameter in the above cited decision tree is, therefore, the determination of the free acidity of the  
53 oil sample; this is defined as the amount of free fatty acids, no longer linked to their parent  
54 triglyceride molecules (TAGs), and measured as percentage of oleic acid. Specifically, the top-  
55 quality product, extra virgin olive oil (EVOO), features a maximum FA of 0.8 g oleic acid/100 g  
56 oil, then the VOO features a maximum FA value of 2.0 g oleic acid/100 g oil and, finally, the  
57 Lampante Olive Oil (LOO), that is not suitable for the commercialization as it is, is characterized by  
58 a FA higher than 2.0 g oleic acid/100 g oil.

59 This parameter is especially affected by the quality of the olives used to produce the oil since free  
60 fatty acids arise from the separation of fatty acids from TAGs because of the action of enzymes,  
61 further stimulated by light, water, and heat. Thus, FA is an indicator of how fresh and how well  
62 handled the olives were before being milled (Tena, Wang, Aparicio-Ruiz, García-González, &  
63 Aparicio, 2015). To confirm this, different studies reported that geographic and environmental  
64 factors (Bustan et al., 2014) and the application of specific technological processes (such as  
65 filtration or a cooling treatment of olive paste) (Veneziani et al., 2018 a-b) do not affect  
66 significantly this parameter.

67 The official procedure to measure the oil free acidity is defined by the European regulations (EEC  
68 Reg. 2568/1991 and following amendments) and consists of an acid-base titration that, albeit simple

77 The possibility of simple, quick and in-situ analysis for the food quality control (often implemented  
78 in the form of portable electronic systems) has been widely researched in the last years (Oujji et al.,  
79 2014; Grossi, Di Lecce, Arru, Gallina Toschi, & Riccò, 2015; Sture, Ruud Oye, Skavhaug, &  
80 Mathiassen, 2016; Arsalane et al., 2018). Concerning the olive oil free acidity determination,  
81 different innovative techniques have been presented in literature (Valli et al., 2016). One approach  
82 is based on Near-Infrared (NIR) spectroscopy, allowing to estimate acidity by optical spectroscopy  
83 analysis of the transmission spectra in the wavenumber range 4541 to 11726  $\text{cm}^{-1}$  (Armenta,  
84 Garrigues, & de la Guardia, 2007). This solution has the advantage that the measurement can be  
85 carried out on the oil sample without any reagent, but it needs expensive instrumentation (optical  
86 spectrophotometer) and requires frequent re-calibration (different calibrations must be carried out  
87 depending on olives varieties or geographical origin etc.). A pH-metric procedure to measure oil  
88 acidity has been also proposed (Tur'yan, Berezin, Kuselman, & Shenhar, 1996): the oil sample is  
89 extracted with a chemical solvent (0.20 M triethanolamine in a 1:1 solution of water and  
90 isopropanol) and the pH measured. This technique is reported to be accurate, but the pH probe  
91 needs to be used with care and frequent cleaning and re-calibration are needed. Moreover, the used  
92 chemicals add costs for disposal after use.

93 The aim of this work is to present a portable electronic system to measure the free acidity in olive  
94 oil samples: the instrument is battery-operated and can be used for quick “in-situ” measurements in  
95 the oil production or bottling site. The working principle is based on the estimation of the oil FA  
96 from the electrical conductance measured by Electrical Impedance Spectroscopy (EIS) of the oil  
97 emulsion with a hydroalcoholic solution. EIS is used in a wide range of applications (Grossi, &  
98 Riccò, 2017a), such as to estimate the ripening degree of fruits (Harker, & Maindonald, 1994), to  
99 characterize plant tissues (Lin, Chen, & Chen, 2012; Ben Hamed, Zorrig, & Hichem Hamzaoui,  
100 2016), to characterize and to detect the freezing end point of ice cream mixes (Grossi, Lanzoni,  
101 Lazzarini, & Riccò, 2012a), to estimate the water content in extra virgin olive oil (Ragni et al.,  
102 2013), to investigate the corrosion of metal surfaces exposed to acid electrolytes (Ansari, Quraishi,

111 instrument has been in-house validated evaluating its performances in terms of correlation between  
112 the concentration of the analyte and the instrument response, limit of detection (LOD) and  
113 quantification (LOQ), precision and accuracy. The validation step, in fact, is one of the measures  
114 universally recognized as a necessary part of a comprehensive system of quality assurance in  
115 analytical chemistry and is an essential component of the measures that a laboratory should  
116 implement to allow it to produce reliable analytical data (Thompson, Ellison, & Wood, 2002).

117

## 118 **2. Material and methods**

### 119 **2.1 Reagents and chemicals**

120 Diethyl ether (ACS reagent, purity  $\geq 99.8\%$ ), ethanol (ACS reagent, purity  $\geq 96\%$ ) and  
121 phenolphthalein solution (indicator, 1% in ethanol) were supplied from Sigma-Aldrich, Inc. (St.  
122 Louis, MO, USA). Sodium hydroxide 0.1 mol/L (N/10) was purchased from Carlo Erba Reagents  
123 S.r.l. (Milan, Italy).

124 Distilled water was produced by Elix Essential system (Millipore, Molsheim, France). Oleic acid  
125 (CAS Number 112-80-1, assay 90%) used as standard for building the calibration curve was  
126 supplied from Sigma-Aldrich, Inc. (St. Louis, MO, USA).

### 127 **2.2 Samples**

128 In order to calibrate the system, 5 samples of refined sunflower oil (with free acidity and peroxide  
129 value under the limits for vegetables oils according to CODEX STAN 210-1999) added with oleic  
130 acid were prepared. To cover the range of free acidity of the different quality grades (EVOO, VOO,  
131 LOO) different concentrations of oleic acid (from 0.25% to 3.75%) were added (Table 1).

132 The in-house validation study was carried out on a set of 30 olive oil (numerically coded from 1 to  
133 30) with a different quality grade: 20 EVOOs, 7 VOOs and 3 LOOs. An aliquot of 50 mL for each  
134 sample was stored in PET dark bottles at 11-12 °C before the analysis.

### 135 **2.3 Determination of free acidity**

#### 136 **2.3.1 Official method**

145 the hydroalcoholic solution, the free fatty acid molecules  $\text{RCOOH}$ , where R is the hydrocarbon  
146 chain, dissociate in the ionic compounds  $\text{H}_3\text{O}^+$  and  $\text{RCO}_2^-$  that contribute to the increase of the  
147 emulsion electrical conductance. Consequently, the higher the free fatty acid molecules  
148 concentration, the higher ions concentration and the higher the electrical conductance.

149 A 50 mL round bottom polypropylene tube (Falcon) modified to feature a couple of cap-shaped  
150 stainless-steel electrodes (6 mm in diameter, spaced by 12 mm one from the other) required for the  
151 electrical characterization is used as sensor to realize the measurement.

152 The first step to realize the measurement is a “reagent test” to check the electrical properties  
153 (electrical conductance) of the hydroalcoholic solution in order to avoid an overestimation of the oil  
154 free acidity linked to the conductance of the solution. For this purpose, the tube (sensor) is filled  
155 with 9 mL of ethanol and 6 mL of distilled water and its conductance is measured. If the reagent  
156 passes the test and is suitable for the measure (conductance not higher than  $0.6 \mu\text{S}$ ), 1 mL of the oil  
157 sample to be tested is added **using a plastic graduated Pasteur pipette** and the obtained mixture is  
158 manually shaken vigorously for about 20 seconds to create an emulsion, whose conductance is  
159 measured to estimate the oil acidity. **In order to avoid error due to the eventual instability of the**  
160 **emulsion, the free acidity measurement should be realized within 30 minutes after the agitation.** A  
161 video showing the procedure to realize the analysis is available in the online version of the paper.

162 The instrument is also equipped with a temperature sensor to measure the room environmental  
163 temperature during the analysis. The emulsion electrical conductance and the temperature are  
164 measured and these values are used to calculate the emulsion conductance at the calibration  
165 temperature ( $23.5 \text{ }^\circ\text{C}$ ). Then the free acidity of the sample is estimated using the calibration  
166 function stored inside the microcontroller non-volatile memory.

#### 167 **2.4 In-house method validation**

168 The following parameters were evaluated in order to ensure the method quality: limit of detection  
169 (LOD) and quantification (LOQ), intra-day and inter-day precision and accuracy.

170 Before the application of the portable system on olive oil samples and in order to verify the



179 concentration of the analyte and the slope of the calibration function ( $m$ ) has been applied. The  
180 formulas used are:

$$181 \text{ LOD} = 3\sigma/m \quad (1)$$

$$182 \text{ LOQ} = 10\sigma/m \quad (2)$$

183 To evaluate the inter-day precision of the method, two virgin olive oil samples for each type of  
184 quality grade (samples EVOO\_2 and EVOO\_19 for EVOOs, samples VOO\_1 and VOO\_6 for  
185 VOOs and samples LOO\_1 and LOO\_3 for LOOs) were analyzed in triplicate for three different  
186 days and values were statistically evaluated by the Student's Test ( $p < 0.05$ ).

187 Accuracy and intra-day precision were studied analyzing the set of 30 olive oils with a different  
188 quality grade. Precision was expressed in terms of RSD (Relative Standard Deviation) measured for  
189 the free acidity values obtained with the portable system. For the accuracy, all the samples were  
190 also analyzed following the official method (acid-base titration) and the comparison between the  
191 values obtained by the two methods was performed (two-tailed paired t-test with  $p < 0.05$ ).

192

### 193 3. Results and discussion

#### 194 3.1 Design of the portable system

195 A portable battery-operated electronic system for in-situ measurements of olive oil free acidity has  
196 been built to allow olive oil quality assurance directly in the oil mills or packaging centers.

197 The dimensions of the instrument are 11 x 15 x 5 cm and some pictures are presented in Figure 2a.

198 The instrument can be powered by USB port or using batteries (3 AAA alkaline batteries 1.5 V).

199 The system is composed of an electronic board designed ad-hoc that performs all the operations to  
200 measure the free acidity, a 2 rows 16 columns LCD screen to output the measure results, four  
201 buttons for user interaction and the sensor previously described (see paragraph 2.3.2). The  
202 electronic board is based on the microcontroller STM32L152RCT6A and its schematic is presented  
203 in Figure 2b. The sensor used for the measurements is realized with electrodes of stainless steel, a

204 material less affected by passivation than other materials. The same sensor could be used many

212 The current drawn by the sensor is converted to a sine-wave voltage  $V_{out}(t)$  by means of a current-  
213 to-voltage converter. Given

$$214 \quad V_{in}(t) = V_{M,in} \times \text{sen}(2\pi ft) \quad (3)$$

215 it is

$$216 \quad V_{out}(t) = V_{M,out} \times \text{sen}(2\pi ft + \varphi) \quad (4)$$

217 where  $V_{M,in}$  and  $V_{M,out}$  are the amplitude of the corresponding signals,  $\varphi$  is the phase difference  
218 between the current through the sample and  $V_{in}(t)$ , while  $f$  is the frequency of the test signal (200  
219 Hz).

220 Both  $V_{in}(t)$  and  $V_{out}(t)$  are acquired by the built-in 12-bits ADC inside the microcontroller using a  
221 sampling frequency of 50 kHz and the sine-wave parameters are calculated using the algorithm  
222 previously presented by Grossi et al. (Grossi, Lanzoni, Lazzarini, & Riccò, 2012b).

223 The emulsion in direct contact with the electrodes can be modeled as the parallel of an electrical  
224 conductance (accounting for the conductance of the emulsion) and a capacitance (accounting for the  
225 emulsion dielectric properties): while the emulsion conductance (that dominates at low frequency)  
226 is affected by the sample acidity due to the variation of the ions concentration, the dielectric  
227 properties are almost independent.

228 The electrical conductance ( $G_m$ ) is thus calculated as:

$$229 \quad G_m = \frac{I}{R_F} \times \frac{V_{M,out}}{V_{M,in}} \times \cos(\varphi) \quad (5)$$

230 where  $R_F$  is the feedback resistance (470 k $\Omega$ ) of the current-to-voltage converter.

231 However, the relation between  $G_m$  and the free acidity is non-linear and it can be modelled with the  
232 function:

$$233 \quad G_m = \alpha + \beta \times \sqrt{FA} \quad (6)$$

234 where  $\alpha$ ,  $\beta$  are empirical parameters that must be determined by a suitable calibration procedure and  
235 are also function of the calibration temperature.

236 Then, the sample free acidity can be estimated from the measured electrical conductance of the

246 emulsion conductance at the calibration temperature ( $G_{m,Tcalib}$ ). Then the free acidity of the oil  
247 sample is estimated from the calculated value of  $G_{m,Tcalib}$  using the calibration function stored inside  
248 the microcontroller non-volatile memory.

249 A more detailed description of the system by an electrical and mathematical point of view is  
250 presented in Appendix A, while a discussion on the influence of the environmental temperature on  
251 the measured conductance and the method used for compensation is presented in Appendix B.

### 252 **3.2 In-house method validation**

253 Before being used with real olive oil samples, the portable instrument has been calibrated. For this  
254 purpose, all the samples of refined sunflower oil added with oleic acid (samples A-E) were analyzed  
255 by both the official method to determine the reference value of free acidity and by the portable  
256 system to measure the conductance value. The data obtained are shown in Table 1. Then the values  
257 were plotted (Figure 3): in all cases the electrical conductance measured for the samples increases  
258 with its free acidity with a non-linear relation, confirming what previously presented (Grossi, Di  
259 Lecce, Gallina Toschi, & Riccò, 2014b). Starting from these data and applying a nonlinear  
260 regression, the equation of the calibration curve was defined as follow:

$$261 \quad FA = \left( \frac{G_{m,23.5^{\circ}C} + 0.0678}{2.7877} \right)^2 \quad (8)$$

262 where  $G_{m,23.5^{\circ}C}$  is the emulsion electrical conductance at the calibration temperature of 23.5 °C.

263 It allows to obtain, directly on the display of the portable system, the free acidity value of the tested  
264 sample estimated starting from its electrical properties. **Since this procedure is necessary to define  
265 the relation between the oil FA and the conductance, and, consequently, the mathematical function  
266 that has to be stored in the microcontroller non-volatile memory of the instrument, this step needs to  
267 be realized only when the instrument is built.**

268 Subsequently, some parameters for the in-house validation of the instrument were evaluated in  
269 order to check its performances and applicability on real olive oil samples.

278 applying this new analytical approach in order to study its accuracy and intra-day precision  
279 (expressed in terms of RSD).

280 The portable instrument showed good precision results since all the obtained values were under the  
281 15% that is considered as acceptable for the validation of a new method (Peters, Drummer, &  
282 Musshoff, 2007). All the samples were also analyzed following the official method (acid-base  
283 titration) and the comparison between the values obtained by the two methods was performed to  
284 measure the accuracy of the system (Table 2). The differences between the two series of results  
285 (official method vs portable system) were evaluated by using the two-tailed paired t-test ( $p < 0.05$ ).  
286 The two approaches did not give statistically differences for the mean values ( $t < t$  (critical value)).  
287 Moreover, the regression between the two series of values provided a coefficient  $R^2_{adj}$  of 0.97  
288 (Figure 4) in agreement with data previously presented by Grossi et al. (Grossi, Di Lecce, Gallina  
289 Toschi, & Riccò, 2014). Considering the commercial categories of the samples analyzed, all of  
290 them, with the only exception of the sample VOO\_7, were classified in the same way by both  
291 approaches. Finally, considering the inter-day precision of the instrument (Table 3), no significative  
292 differences (Student's Test,  $p < 0.05$ ) were found among the results obtained for each analyzed  
293 sample in the three different days.

294

#### 295 4. Conclusions

296 The design and in-house validation of a portable battery-operated electronic system suitable for in-  
297 situ measurements of olive oil free acidity has been presented. The system is built with low cost  
298 electronics and embeds a temperature sensor to compensate variations of the measured electrical  
299 parameters with the environmental temperature, thus making it suitable for on-site free acidity  
300 measurements outside a laboratory. Its working principle is based on the estimation of the olive oil  
301 FA from the measure of the conductance of an emulsion between a hydro-alcoholic solution and the  
302 sample to be tested. When the free fatty acids present in the sample come in contact with the  
303 hydroalcoholic solution, a dissociation occurs, leading to the formation of ions that produce an

312 focused on the evaluation of the performance of this portable system when applied to bottles stored  
313 for several months in dark/light conditions simulating possible different commercial  
314 conservation/exposure on shelves of groceries. This study will give indications regarding the  
315 maximum oxidation state of olive oil sample tolerable to use this portable system as accurate and  
316 reliable.

317

### 318 **Funding**

319 This work was supported by the Horizon 2020 European Research project OLEUM “*Advanced*  
320 *solutions for assuring the authenticity and quality of olive oil at a global scale*”, which has received  
321 funding from the European Commission within the Horizon 2020 Programme (2014–2020), grant  
322 agreement no. 635690. The information expressed in this article reflects the authors’ views; the  
323 European Commission is not liable for the information contained herein.

324

### 325 **Declaration of interest**

326 None.

327

### 328 **References**

329

330 Ansari, K. R., Quraishi, M. A., Singh, A. (2015). Pyridine derivatives as corrosion inhibitors for  
331 N80 steel in 15% HCl: electrochemical, surface and quantum chemical studies. *Measurement*, 76,  
332 136-147.

333

334 Armenta, S., Garrigues, S., de la Guardia, M. (2007). Determination of edible oil parameters by  
335 near infrared spectrometry. *Analytica Chimica Acta*, 596, 330-337.

336

337 Arsalane A, El Barbri N, Tabyaoui A, Klilou A, Rhofir K, Halimi A (2018). An embedded

345 Bendini, A., Cerretani, L., Carrasco-Pancorbo, A., Gomez-Caravaca, A. M., Segura-Carretero, A.,  
346 Fernandez-Gutierrez, A., Lercker, G. (2007). Phenolic molecules in virgin olive oils: a survey of  
347 their sensory properties, health effects, antioxidant activity and analytical methods. An overview of  
348 the last decade. *Molecules*, *12*, 1679-1719.

349

350 Bera, T. K. (2014). Bioelectrical impedance methods for noninvasive health monitoring: a review.  
351 *Journal of Medical Engineering*, article ID 381251, 1-28.

352

353 Bustan, A., Kerem, Z., Yermiyahu, U., Ben-Gal, A., Lichter, A., Droby, S., Zchori-Fein, E.,  
354 Orbach, D., Zipori, I., Dag, A. (2014). Preharvest circumstances leading to elevated oil acidity in  
355 'Barnea' olives. *Scientia Horticulturae*, *176*, 11-21.

356

357 Codex Alimentarius – CODEX STAN 210-1999 - Standard for named vegetable oils codex.  
358 Adopted in 1999. Revised in 2001-2003-2009-2017. Amended in 2005-2011-2013-2015.

359

360 Commission Implementing Regulation (EU) No 2016/1227 of 27 July 2016 amending Regulation  
361 (EEC) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the relevant  
362 methods of analysis.

363

364 Commission Implementing Regulation (EU) No 1348/2013 of 16 December 2013 amending  
365 Regulation (EEC) No 2568/91 on the characteristics of olive oil and olive-residue oil and on the  
366 relevant methods of analysis.

367

368 Commission Regulation (EEC) No 2568/91 of 11 July 1991 on the characteristics of olive oil and  
369 on the relevant methods of analysis.

370

378 Grossi, M., Lanzoni, M., Lazzarini, R., Riccò, B. (2012b). Linear non-iterative sinusoidal fitting  
379 algorithm for microbial impedance biosensor. *Sensors & Transducers Journal*, 137 (2), 235-244.  
380

381 Grossi, M., Di Lecce, G., Gallina Toschi, T., Riccò, B. (2014a). A novel electrochemical method  
382 for olive oil acidity determination. *Microelectronics Journal*, 45 (12), 1701-1707.  
383

384 Grossi, M., Di Lecce, G., Gallina Toschi, T., Riccò, B. (2014b). Fast and accurate determination of  
385 olive oil acidity by electrochemical impedance spectroscopy. *IEEE Sensors Journal*, 14 (9), 2947-  
386 2954.  
387

388 Grossi, M., Di Lecce, G., Arru, M., Gallina Toschi, T., Riccò, B. (2015). An opto-electronic system  
389 for in-situ determination of peroxide value and total phenol content in olive oil. *Journal of Food*  
390 *Engineering*, 146, 1-7.  
391

392 Grossi, M., Riccò, B. (2017a). Electrical impedance spectroscopy (EIS) for biological analysis and  
393 food characterization: a review. *Journal of Sensors and Sensor Systems*, 6, 303-325.  
394

395 Grossi, M., Riccò, B. (2017b). An automatic titration system for oil concentration measurement in  
396 metalworking fluids. *Measurement*, 97, 8-14.  
397

398 Harker, F. R., Maindonald, J. H. (1994). Ripening of nectarine fruit. *Plant Physiology*, 106, 165-  
399 171.  
400

401 Hussin, M. H., Rahim, A. A., Nasir, M., Ibrahim, M., Brosse, N. (2016). The capability of  
402 ultrafiltrated alkaline and organosolv oil palm (*Elaeis guineensis*) fronds lignin as green corrosion  
403 inhibitor for mild steel in 0.5 M HCl solution. *Measurement*, 78, 90-103.

412 Lin, C. M., Chen, L. H., Chen, T. M. (2012). The development and application of an electrical  
413 impedance spectroscopy measurement system for plant tissue. *Computers and Electronics in*  
414 *Agriculture*, 82, 96-99.

415

416 Oujji, N. B., Bakas, I., Istamboulié, G., Ait-Ichou, I., Ait-Addi, E., Rouillon, R., Noguier, T. (2014).  
417 A simple colorimetric enzymatic-assay, based on immobilization of acetylcholinesterase by  
418 adsorption, for sensitive detection of organophosphorus insecticides in olive oil. *Food Control*, 46,  
419 75-80.

420

421 Peters, F. T., Drummer, O. H., Musshoff F. (2007). Validation of new methods. *Forensic Science*  
422 *International*, 165, 216–224.

423

424 Ragni, L., Iaccheri, E., Cevoli, C., Berardinelli, A., Bendini, A., Gallina Toschi, T. (2013). A  
425 capacitive technique to assess water content in extra virgin olive oils. *Journal of Food Engineering*,  
426 116 (1), 246-252.

427

428 Sture, O., Ruud Oye, E., Skavhaug, A., Mathiassen, J. R. (2016). A 3D machine vision system for  
429 quality grading of Atlantic salmon. *Computers and Electronics in Agriculture*, 123, 142-148.

430

431 Tena, N., Wang, S. C., Aparicio-Ruiz, R., García-González, D. L., Aparicio R. (2015). In-depth  
432 assessment of analytical methods for olive oil purity, safety, and quality characterization. *Journal of*  
433 *Agricultural and Food Chemistry*, 63, 4509–4526.

434


435 Thompson, M., Ellison, S. L., Wood, R. (2002). Harmonized guidelines for single-laboratory  
436 validation of methods of analysis (IUPAC Technical Report). *Pure and Applied Chemistry*, 74 (5),  
437 835-855.




447 Veneziani, G., Esposito, S., Taticchi, A., Urbani, S., Selvaggini, R., Sordini, B., Servili, M. (2018a).  
448 Characterization of phenolic and volatile composition of extra virgin olive oil extracted from six  
449 Italian cultivars using a cooling treatment of olive paste. *LWT – Food Science and Technology*, 87,  
450 523-528.

451


452 Veneziani, G., Esposito, S., Minnocci, A., Taticchi, A., Urbani, S., Selvaggini, R., Sordini, B.,  
453 Sebastiani, L., Servili, M. (2018b). Compositional differences between veiled and filtered virgin  
454 olive oils during a simulated shelf life. *LWT – Food Science and Technology*, 94, 87-95.

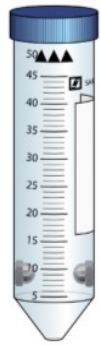

The sensor is filled with the hydroalcoholic solution (40% distilled water / 60% ethanol) and, after the addition of 1 ml of the oil sample, is manually shaken to create the emulsion

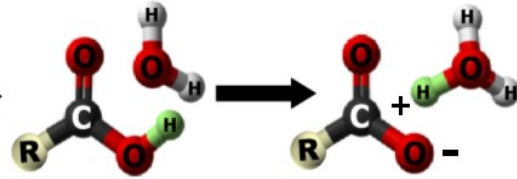

The free fatty acid molecules dissociate in presence of the hydroalcoholic solution and the  $\text{H}_3\text{O}^+$  ions contribute to the increase of the emulsion electrical conductance


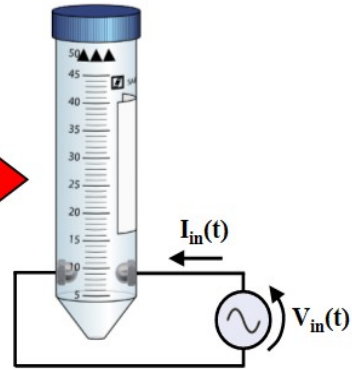


A sine-wave signal (1V amplitude, 200 Hz frequency) is applied to the sensor electrodes. The emulsion electrical conductance is measured and the oil sample free acidity is calculated

455  
456

Figure 1

(a)

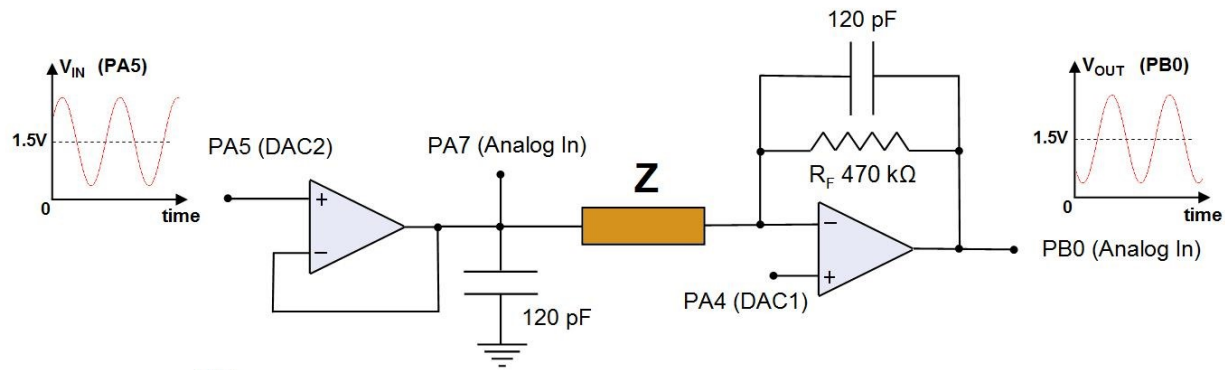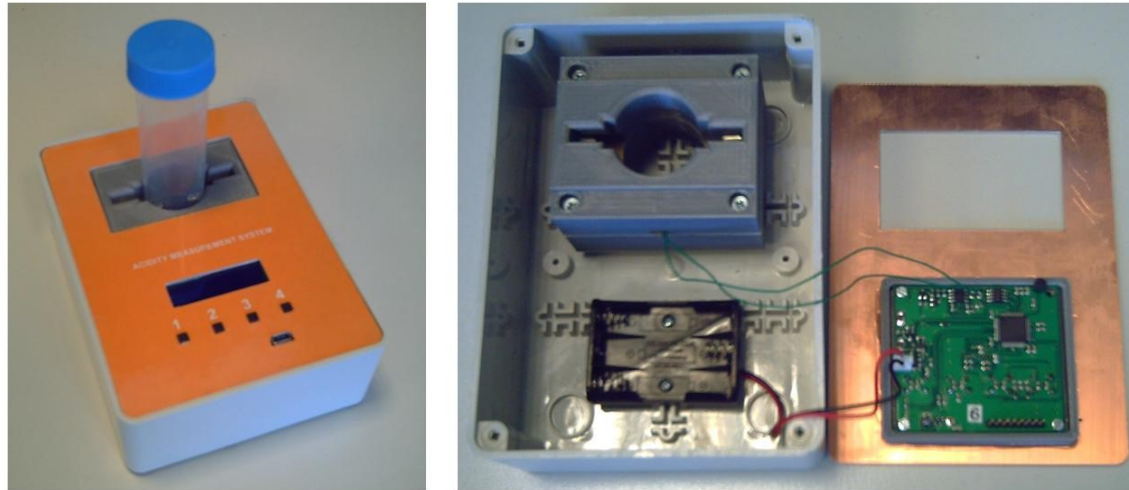

(b)


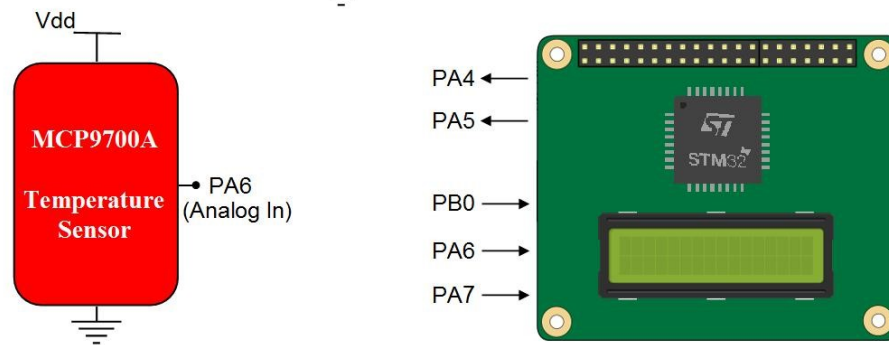


Figure 2

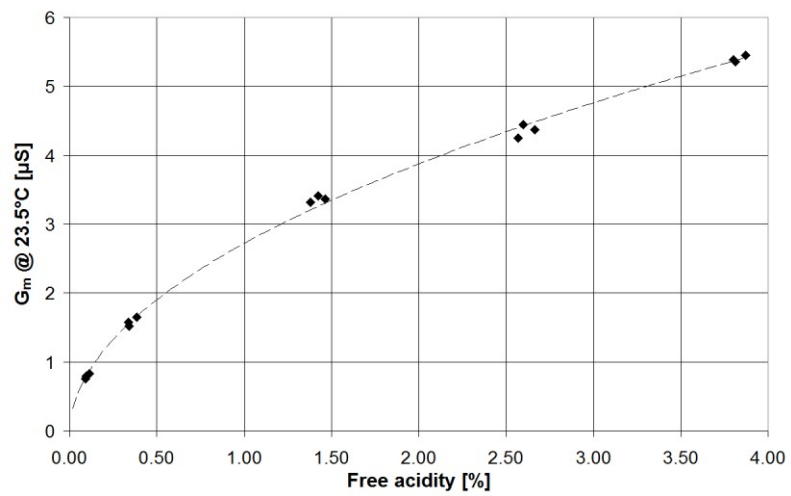


460

461

462

**Figure 3**

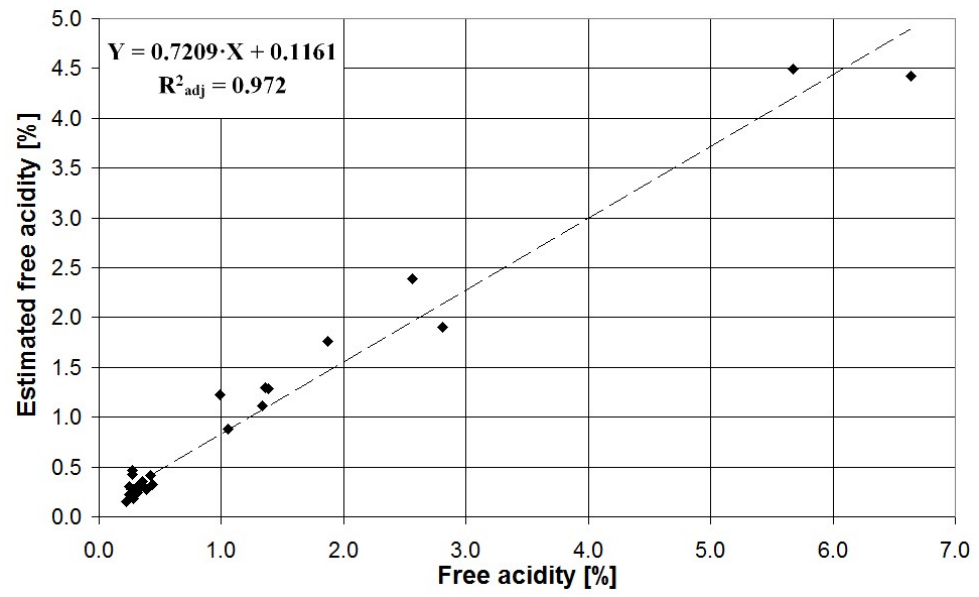

463

464

Figure 4

465 **Figure captions**

466 **Figure 1.** Scheme of the working principle of the proposed electronic system for the determination  
467 of free acidity in virgin olive oil.

468 **Figure 2.** Pictures of the outside and inside of the electronic system and the sensor for olive oil free  
469 acidity analysis (a); system hardware and electrical scheme (b).

470 **Figure 3.** Measured electrical conductance ( $G_m$ ) at 23.5°C vs free acidity for the sunflower oil  
471 calibration set.

472 **Figure 4.** Scatter plot of the estimated free acidity vs the free acidity measured by titration for a set  
473 of 30 olive oil samples.

474 **Table 1**

475 Codes and description of the samples used for the calibration of the system, data of free acidity  
476 determined by the official method and conductance measured with the portable system.

<b>Sample code</b>	<b>Sample description</b>	<b>Free acidity (% oleic acid)</b>	<b>Conductance (<math>\mu</math>S)</b>
A	Refined sunflower oil (without addition of oleic acid)	$0.10 \pm 0.01$	$0.79 \pm 0.03$
B	Refined sunflower oil + 0.25% oleic acid	$0.36 \pm 0.03$	$1.58 \pm 0.06$
C	Refined sunflower oil + 1.25% oleic acid	$1.4 \pm 0.04$	$3.36 \pm 0.04$
D	Refined sunflower oil + 2.50% oleic acid	$2.6 \pm 0.05$	$4.35 \pm 0.10$
E	Refined sunflower oil + 3.75% oleic acid	$3.8 \pm 0.04$	$5.40 \pm 0.05$

477 The values are mean of three replicates. According to the EU Reg. 2016/1227, the FA values from 0  
478 up to 1 (including 1) are reported with two decimal places and the FA values higher than 1 are  
479 reported with one decimal place.

480 **Table 2.** Values of free acidity (mean and standard deviation) for all samples measured by the  
 481 portable system and the official method.

Sample code	Free acidity (portable system)	Free acidity (official method)
EVOO_1	0.18 ± 0.01	0.25 ± 0.01
EVOO_2	0.31 ± 0.03	0.34 ± 0.02
EVOO_3	0.25 ± 0.01	0.28 ± 0.01
EVOO_4	0.43 ± 0.01	0.27 ± 0.01
EVOO_5	0.30 ± 0.02	0.25 ± 0.01
EVOO_6	0.22 ± 0.01	0.25 ± 0.01
EVOO_7	0.18 ± 0.02	0.25 ± 0.01
EVOO_8	0.29 ± 0.04	0.37 ± 0.02
EVOO_9	0.29 ± 0.01	0.34 ± 0.01
EVOO_10	0.47 ± 0.01	0.27 ± 0.01
EVOO_11	0.32 ± 0.02	0.33 ± 0.01
EVOO_12	0.41 ± 0.02	0.42 ± 0.02
EVOO_13	0.24 ± 0.02	0.28 ± 0.01
EVOO_14	0.27 ± 0.03	0.28 ± 0.01
EVOO_15	0.32 ± 0.02	0.44 ± 0.02
EVOO_16	0.35 ± 0.02	0.35 ± 0.01
EVOO_17	0.18 ± 0.02	0.28 ± 0.01
EVOO_18	0.15 ± 0.01	0.22 ± 0.01
EVOO_19	0.27 ± 0.02	0.39 ± 0.02
EVOO_20	0.24 ± 0.00	0.31 ± 0.01
VOO_1	1.3 ± 0.09	1.4 ± 0.02
VOO_2	1.3 ± 0.02	1.4 ± 0.03
VOO_3	0.88 ± 0.07	1.1 ± 0.01
VOO_4	1.1 ± 0.05	1.3 ± 0.00
VOO_5	1.2 ± 0.05	0.99 ± 0.00
VOO_6	1.8 ± 0.12	1.9 ± 0.01
VOO_7	1.9 ± 0.07	2.8 ± 0.02
LOO_1	2.4 ± 0.32	2.6 ± 0.00
LOO_2	4.4 ± 0.04	6.6 ± 0.05
LOO_3	4.5 ± 0.17	5.7 ± 0.01

482 Results are expressed as % of oleic acid. According to the EU Reg. 2016/1227, the FA values from  
 483 0 up to 1 (including 1) are reported with two decimal places and the FA values higher than 1 are



485 **Table 3.** Values of free acidity (mean and standard deviation) measured in three different days to  
486 evaluate the inter-day precision.

<b>Sample</b>	<b>Day 1</b>	<b>Day 2</b>	<b>Day 3</b>
EVOO_2	0.31 ± 0.03	0.34 ± 0.01	0.32 ± 0.02
EVOO_19	0.27 ± 0.02	0.25 ± 0.01	0.27 ± 0.03
VOO_1	1.3 ± 0.09	1.3 ± 0.02	1.3 ± 0.01
VOO_6	1.8 ± 0.12	1.8 ± 0.08	1.8 ± 0.05
LOO_1	2.4 ± 0.32	2.5 ± 0.18	2.5 ± 0.11
LOO_3	4.5 ± 0.17	4.6 ± 0.25	4.4 ± 0.04

487 EVOO\_2 and EVOO\_19: extra virgin olive oils; VOO\_1 and VOO\_6: virgin olive oils; LOO\_1 and  
488 LOO\_3: lampante olive oils. According to the EU Reg. 2016/1227, the FA values from 0 up to 1  
489 (including 1) are reported with two decimal places and the FA values higher than 1 are reported  
490 with one decimal place.

491 **Appendix A**

492 According to the working principle of the portable system, a 50 mL polypropylene tube (Falcon)  
493 vial modified to feature a couple of stainless-steel electrodes to measure the emulsion conductance  
494 (hereafter the sensor) is filled with 15 mL of hydro-alcoholic solution (40% distilled water/60%  
495 ethanol), then 1 mL of the olive oil sample is added and all is stirred to create an emulsion.

496 In presence of the hydroalcoholic solution, the free fatty acid molecule RCOOH, where R is the  
497 hydrocarbon chain, dissociates in the ionic compounds  $\text{H}_3\text{O}^+$  and  $\text{RCO}_2^-$  that contribute to the  
498 increase of the emulsion electrical conductance. In the end, the higher the free fatty acid molecules  
499 concentration, the higher ions concentration and the higher the electrical conductance.

500 The emulsion electrical conductance is measured by Electrical Impedance Spectroscopy (EIS).

501 In the proposed approach the sample under investigation is stimulated with a sine-wave voltage  
502 signal  $V_{in}(t)$ :

503 
$$V_{in}(t) = V_{M,in} \times \sin(2\pi ft) \quad (\text{A1})$$

504 and the current  $I_{in}(t)$  through the sample is measured:

505 
$$I_{in}(t) = I_{M,in} \times \sin(2\pi ft + \varphi) \quad (\text{A2})$$

506 where  $V_{M,in}$  and  $I_{M,in}$  are the amplitudes of the corresponding signals,  $f$  is the frequency of the test  
507 signal and  $\varphi$  is the phase difference between  $I_{in}(t)$  and  $V_{in}(t)$ .

508 The sample electrical admittance is then expressed as:

509 
$$Y = \frac{I_{in}(j2\pi f)}{V_{in}(j2\pi f)} = \frac{I_{M,in}}{V_{M,in}} \times (\cos \varphi + j \times \sin \varphi) = \text{Re}(Y) + j \times \text{Im}(Y) \quad (\text{A3})$$


510 The emulsion in direct contact with the electrodes can be modeled as the parallel of an electrical  
511 conductance (accounting for the conductance of the emulsion) and a capacitance (accounting for the  
512 emulsion dielectric properties): while the emulsion conductance (that dominates at low frequency)  
513 is affected by the sample acidity due to the variation of the ions concentration, the dielectric  
514 properties are almost independent. Thus, the emulsion electrical conductance  $G_m$  can be estimated  
515 with the real component of the emulsion admittance  $\text{Re}(Y) = |Y| \cdot \cos(\varphi)$  where  $|Y|$  is the admittance

524 
$$G_m = \alpha + \beta \times \sqrt{FA} + \gamma \times \sqrt[4]{FA^3} \quad (A4)$$

525 where  $\alpha$ ,  $\beta$  and  $\gamma$  are empirical parameters that must be determined by a suitable calibration  
 526 procedure and are also function of the calibration temperature.

527 In Figure A1 (b) the qualitative plot of  $G_m$  vs FA is shown: as can be seen, the non-linear function  
 528 results in better accuracy for the estimated free acidity for lower acidity levels. This has been taken  
 529 in account by using a suitable amount of sample to create the emulsion to obtain good accuracy in  
 530 the acidity range of interest.

531

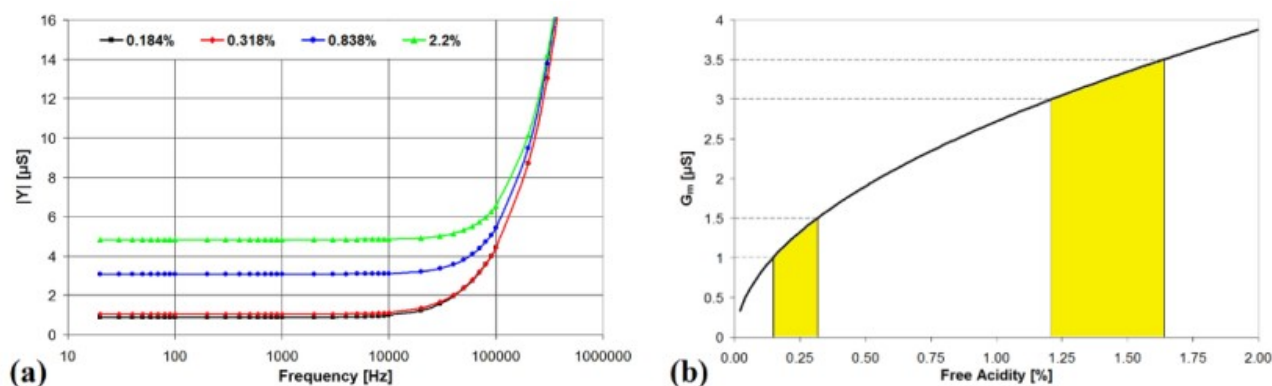

532 **Figure A1** (a) admittance modulus plotted vs frequency for olive oil samples featuring different free  
 533 acidity; (b) qualitative plot of the electrical conductance as function of sample acidity.

534 (permission for publishing the Figure A1(a) obtained from Grossi et al., 2014. *Microelectronics*  
 535 *Journal*, 45 (12), 1701-1707)

536

537  
 538 However, eq. A4 needs a computation intensive iterative algorithm to extract the estimated acidity  
 539 from the measured value of  $G_m$  as well as manual input of starting point to avoid failing in  
 540 algorithm convergence. Thus, a simpler model has been chosen (obtained by neglecting the molar  
 541 conductivity dependence on the  $H_3O^+$  ions concentration) that is much more suitable to be  
 542 implemented in a low-cost microcontroller and can be described by the following function:

550


## 551 **Appendix B**

552 The oil sample free acidity can be estimated by measuring the emulsion electrical conductance at  
553 the temperature of calibration and then calculating the free acidity using equation A6. However,  
554 since the system must be operated “in the field” and the environmental temperature is not a  
555 parameter under control, there is the need to investigate how the emulsion conductance varies with  
556 the temperature so that the oil free acidity can be estimated by the measure of the emulsion  
557 conductance and the temperature.

558 Four different olive oil samples featuring different free acidity values (0.06% sample A, 0.37%  
559 sample B, 0.81% sample C and 3.7% sample D) were tested inside a Binder APT KB 53 thermal  
560 incubator for different temperatures between 15°C and 35°C.

561 In Figure B1 (a) the measured emulsion conductance is plotted vs the incubation temperature for  
562 each sample. In all cases the  $G_m$  is a linear function of the temperature with determination  
563 coefficients  $R^2 > 0.99$ . The calculated linear regression lines allow to determine the conductance  
564 variation with temperature (i.e.  $\partial G_m / \partial T$ ) for all samples: 0.0196 for sample A, 0.0495 for sample B,  
565 0.0579 for sample C and 0.1417 for sample D.  $\partial G_m / \partial T$  is thus found to increase with the sample  
566 free acidity. Since the sample free acidity is also a function of the emulsion electrical conductance,  
567 the relation between  $\partial G_m / \partial T$  and  $G_{m,23.5^\circ\text{C}}$  has been plotted in Figure B1 (b).

568

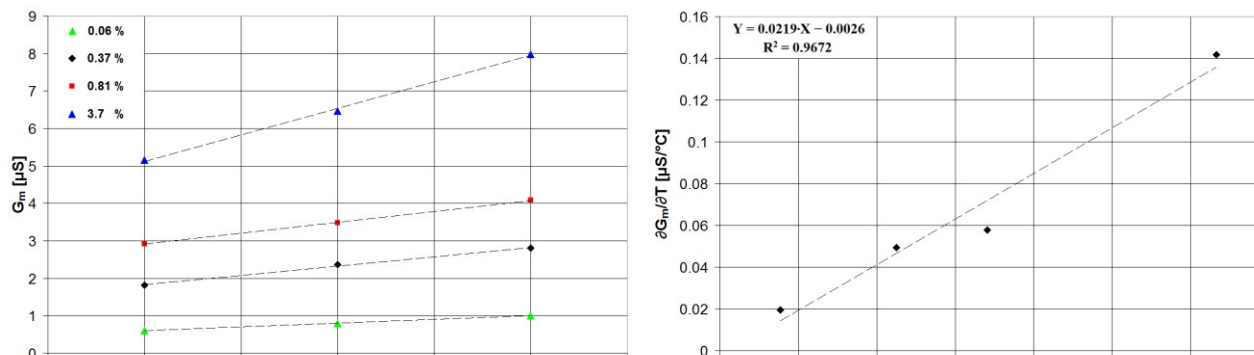

576 where  $\partial G_m/\partial T$  and  $G_{m,23.5^\circ C}$  are expressed as  $\mu S/^\circ C$  and  $\mu S$  respectively. The emulsion electrical  
577 conductance at the environmental temperature  $T$  can thus be expressed as:

$$578 \quad G_{m,T} = G_{m,23.5^\circ C} + \frac{\partial G_m}{\partial T} \times (T-23.5) \quad (B2)$$

579 and

$$580 \quad G_{m,T} = G_{m,23.5^\circ C} + (0.0219 \times G_{m,23.5^\circ C} - 0.0026) \times (T-23.5) \quad (B3)$$

581 The emulsion electrical conductance at  $T_{calib} = 23.5^\circ C$  can thus be estimated from the electrical  
582 conductance at temperature  $T$  and the measured value of  $T$  using the following formula:

$$583 \quad G_{m,23.5^\circ C} = \frac{G_{m,T} + 0.0026 \times (T-23.5)}{1 + 0.0219 \times (T-23.5)} \quad (B4)$$

584 Thus, by measuring  $G_{m,T}$  and  $T$ , the value of  $G_{m,23.5^\circ C}$  can be calculated using equation B4 and, from  
585 this value, the sample acidity can be estimated using equation A6.